

SCENAR COSMODIC EX715Ag

Personal Devices



SCENAR COSMODIC EX715Ag is an exclusive SCENAR-COSMODIC series model employing COSMODIC technology. Intended for personal use. Completely automated treatment process and basic parameters.

Innovations

This is the first EX model which was endowed with OLED display. Its sophisticated design allows anyone operate the device easily. Even battery replacement is like opening a jewel-box — all you need is press the little button.

Regulation of action energy level is now performed automatically.

SCENAR COSMODIC EX715Ag automatically switches to the more delicate impulse for action upon mucous membranes (when you work with vaginal and rectal electrodes).

The graphics in area search mode are simple to follow and easy to understand.

The menu allows you to set the user's functions, in Russian and English

Only with this model you can enjoy the benefits of 3 cosmetological modes.

The quiet mode is implemented for the first time: you switch-off sound and activate vibrosignal to indicate the reaction. Now you can treat yourself during a staff meeting or in a theatre without disturbing anyone (a toothache or headache may bother you anywhere).

Additional option is allowed for — electronic reference book.

This device is the first one which implements COSMODIC technology in full.

Main features

- “dynamic adaptation” block allows you to automatically manage stimulation during prolonged action
- four action modes: one universal mode for first aid and treatment and three modes for cosmetic purposes to act upon face, middle part of the body and legs
- on-screen menu allows you to manage the duration of cosmetic procedures, switch off the sound and set the language (Russian/English)
- interactive help with pictures allows you to quickly find appropriate areas for action
- operating system lets you simultaneously look through the interactive help while you are conducting treatment, without interrupting the action
- updating interactive help via Internet
- timer of action ,procedure, & resource of the device.
- graphic gray-shade OLED display
- LED and sound indication;
- vibrosignal;
- detachable silver electrodes;
- connecting remote electrode;.

Touch of Power Therapies

PO Box 239

Mapleton

Queensland 4560

Australia

www.touchofpower.com.au

ph 07 5445 7778

m 0407 136 794

Indication

Graphic gray-shade OLED display

Shows action modes, action energy, activity of the area being acted upon, dosage, time of action upon the certain area and full time of action, menu, interactive help, and battery charge level.


LED and sound indication

Allows you to approximately define appropriate areas for action, activation of dynamic adaptation block and proper contact of the electrodes with skin.

Interactive help is a short guide with pictures which allow you to quickly find proper areas for action. Interactive help is an essential thing when you have no detailed instructions at hand, since it is not always convenient to carry the book of instructions with you. The operating system allows you to use the interactive help, look through the pictures and come back to the procedure displays even during the treatment without removing the electrodes from the skin surface. The displays of treatment and the picture you are looking at in the help are saved during the procedure, so that you could quickly shift from the procedure displays to the picture. The help may consist of several directories. You will see a list of directories on the screen (the number of directories depends on your order). Below is an example of a directory you may see on the screen of your device.

Directory “Systems” consists of three levels:

- list of nosological forms,
- list of diseases of each nosological form,
- pictures showing areas which need to be treated.

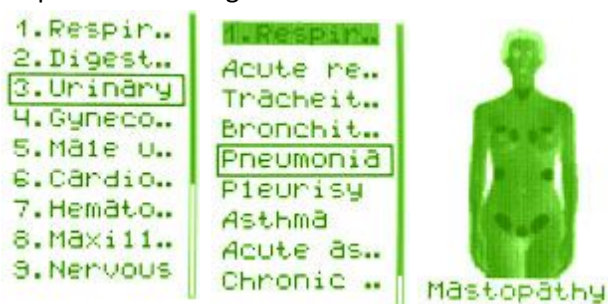

Pictures appear in the order the areas need to be treated. Interactive help is recorded on the RS-MMC flash memory. It allows updating it via Internet.

Electrodes

The electrodes are detachable and made in a form of three bars, which are flat on the under side and rounded at the front end of the device. The rounded part is convenient for treatment of small or hard-to-reach areas (for example, on the face). The flat part is mainly used for treatment of large areas. The electrodes are made of silver according to the special technology, which allows maintaining the natural structure of silver. This provides additional healing effect and more accurate signal of action. Silver is known for its bactericidal, antiseptic, anti-inflammatory and styptic effects. Silver is a natural bactericidal metal effective against 650 kinds of bacteria which do not acquire resistance to it, like to most antibiotics.

In addition, you can order electrodes made of special alloy of gold. Gold electrodes have the same shape as silver ones and can be easily installed in the device by the user.

Casing

The device is encased in a handy and ergonomic “Smart” casing. Waterproof keyboard and sealants allow easy disinfection and cleaning of the device. The casing is made of shockproof ABS plastic.

Size - 104 x 47 x 28 mm

Power supply - Two rechargeable batteries or galvanic cells, size AAA, 1.5V.

Delivery set - 1. Device 2. User's Guide 3. Package!

Functional and technical abundance enclosed in a compact and ergonomic casing make this device one of the best for personal use. If you wish to broaden the limits for your health and well-being, a top device for yourself, this is what you need. Treat yourself beautifully